

---

# *IDM*

---

## *Instruments*

### **Gelbo Flex Tester**

**Model No:** G0002

**Details:** This apparatus is used to determine whether certain laminations of plastic withstand repetitive strain.

By attaching sample pieces of laminated plastic to the two circular clamping disks, via hose clamps, the specimen is twisted and turned to the specific movements. This process will show a visual result of pinholing and delamination.

### **Applications:**

- Laminated plastic films.

### **Features:**

- 45 rpm
- Stroke : 6"
- Twist : 400° Short stroke
- Twist : 440° Long stroke
- Counter for oscillations
- Totaliser

### **Options:**

- Not applicable

### **Standards:**


- ASTM F 392

### **Connections:**

- Electrical: 220/240 VAC @ 50 HZ

### **Dimensions:**

- **H:** 300mm • **W:** 1100mm • **D:** 350mm
- **Weight:** 45kg

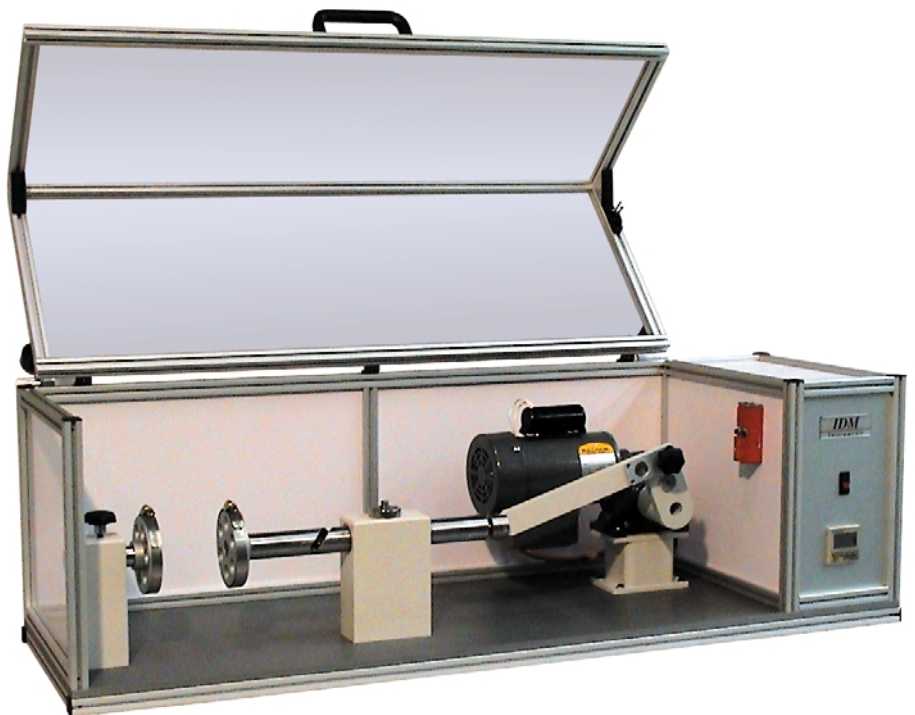

**Representante Oficial para América Latina**  
**REGMED INDÚSTRIA TÉCNICA DE PRECISÃO LTDA.**

tel. +55 11 3601-5700 fax. +55 11 3601-5973 [regmed@regmed.com.br](mailto:regmed@regmed.com.br) [www.regmed.com.br](http://www.regmed.com.br)

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows.  
Products and details shown may change without notice as IDM Instruments reserves the right to make alterations to its product range from time to time.