

---

# *IDM*

---

## *Instruments*

### **Tetrapod Abrasion Tester**

**Model No:** T0004

**Details:** The Tetrapod tester tests the mechanical action on surface texture of textile carpet.

A cylindrical drum, lined with specimens and containing the tetrapod, is rotated the required number of times so that the tetrapod 'walks' over the specimens. The appearance of the tested specimen is compared to that of a concurrently tested controlled specimen.

### **Applications:**

- All textile floor coverings with a thickness no greater than 20mm.

### **Features:**

- 5 digit counter
- Tetrapod: 950 ±20g
- Speed : 50 ±2 rpm
- Horizontal axis of rotation
- Perspex lids

### **Options:**

- Thickness gauge

### **Standards:**


- ISO/TR 6131

### **Connections:**

- Electrical: 220/240 VAC @ 50 HZ  
110 VAC @ 60 HZ

### **Dimensions:**

- **H:** 300mm • **W:** 750mm • **D:** 300mm
- **Weight:** 20kg

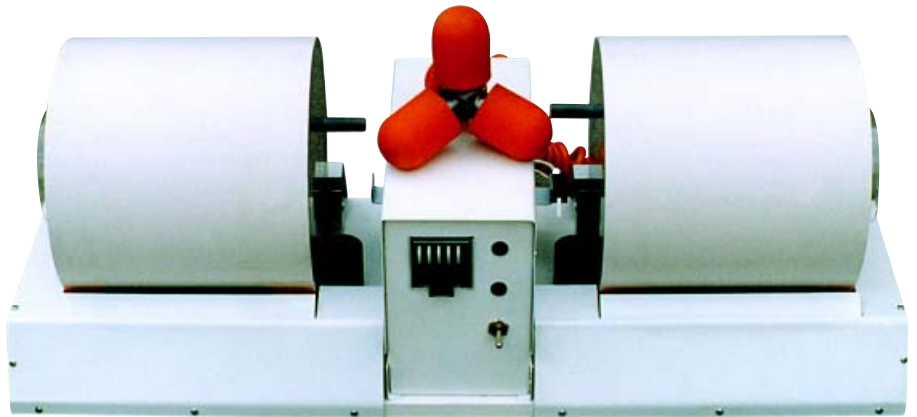

**Representante Oficial para América Latina**  
**REGMED INDÚSTRIA TÉCNICA DE PRECISÃO LTDA.**

tel. + 55 11 3601-5700 fax. + 55 11 3601-5973 [regmed@regmed.com.br](mailto:regmed@regmed.com.br) [www.regmed.com.br](http://www.regmed.com.br)

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows.

Products and details shown may change without notice as IDM Instruments reserves the right to make alterations to its product range from time to time.