

LABORATORY SHEET DRYER TYPE RAPID KÖTHEN model TD-RK



The TD-RK is a laboratory sheet dryer designed to be used along with any Rapid Köthen type Laboratory Sheet Former.

This newly designed dryer meets all specification of applicable standards and replaces traditional Rapid-Köthen drying stations saving laboratory space and maintenance costs.

Its elegant design integrates a balanced combination of simplicity, ruggedness and high technology that ensures excellent performance with minimal requirements for installation, operation and maintenance.

MAJOR TECHNICAL CHARACTERISTICS

Drying area	315 cm ² / Ø 200 mm
Drying temperature - adjustable	(90 to 94)°C
Drying time - adjustable	(0 to 99) min


Rugged construction with simple design uses corrosion resistant materials such as stainless steel, bronze and aluminum. Simple and efficient digital controllers for temperature and drying time.

Optional accessories:

· Standard board carriers - Ø 240mm, 250g/m²

Applicable standards:

· ISO 5269/2

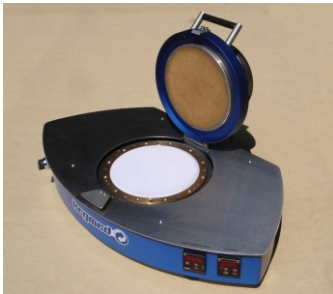

Power supply	220 VAC, single phase, 1200 W
Air supply	2 bar, dry, clean, 0,1 m ³ /h
Dimensions	(65 x 40 x 27) cm (W x L x H)
Weight	45 kg

Note: due to constant development our equipment design and specifications are subject to change without notice.